# KORONA<sup>®</sup>

## GESINA

- D** PERSONENWAAGE Gebrauchsanleitung  
**GB** BATHROOM SCALE Instruction for use  
**FR** PÈSE-PERSONNE Mode d'emploi  
**ES** BÁSCULA DE BAÑO Instrucciones para el uso  
**IT** BILANCIA PESAPERSONE Istruzioni per l'uso  
**GR** ΟΔΗΓΙΕΣ ΧΡΗΣΗΣ Οδηγίες χρήσεως  
**RUS** НАПОЛЬНЫЕ ВЕСЫ Инструкция по применению  
**PL** WAGA OSOBOWA Instrukcja obsługi  
**NL** PERSONENWEEGSCHAAL Gebruikshandleiding  
**PT** BALANÇA DE CASA-DE-BANHO Instruções de utilização  
**TR** TERAZI Kullanma Talimatı  
**DK** PERSONVÆGT Brugsanvisning  
**SV** PERSONVÅG Bruksanvisning  
**FIN** HENKILÖVAAKA Käyttöohje  
**CZ** OSOBNÍ VÁHA Návod k použití  
**SLO** OSEBNA TEHTNICA Navodila za uporabo  
**HU** SZEMÉLYMÉRLEG Használati útmutató

KORONA electric GmbH  
Hauptstraße 169  
59846 Sundern  
Tel.-Nr.: +49 2933 90284-80  
[www.korona-electric.de](http://www.korona-electric.de)

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- GB** Instructions for use
- F** Instructions d'utilisation
- E** Instrucciones de uso
- I** Istruzioni per l'uso
- GR** Οδηγίες χρήσης
- RUS** Инструкция по применению
- PL** Instrukcja obsługi
- NL** Gebruikershandleiding
- P** Instruções de utilização
- TR** Kullanım kılavuzu
- DK** Brugsanvisning
- S** Bruksanvisning
- FIN** Käyttöohje
- CZ** Návod k použití
- SLO** Navodila za uporabo
- H** Használati utasítás

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
## GB General advice

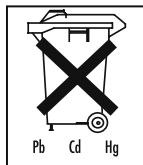- Max. capacity 200 kg (441 lb / 31 st), Graduation 100 g (0.2 lb / 1 lb).
- Cleaning: You can clean the scale with a damp cloth and a little washing up liquid, if required. Never immerse the scale in water or rinse it under running water.
- Do not use aggressive or abrasive cleaning products.
- Protect the scale against knocks, damp, dust, chemicals, drastic changes in temperature, electromagnetic fields and sources of heat that are too near.
- Never step onto the outer edge of the scale with one foot as this could cause the scale to tip!
- Keep packaging material away from children!
- The device is for private use only and is not intended for medical or commercial purposes.
- Should you have any questions regarding operation and usage of our appliances, please contact your local retailers or Customer Service.
- Before making a claim, please check the batteries and replace them if necessary.
- Repairs must only be carried out by Customer Services or authorised suppliers.
- Spent batteries and rechargeable batteries do not constitute normal household waste! They are considered to be toxic waste and, as such, should be disposed of in special containers, toxic waste collection points or brought to electrical goods dealers. You are legally obliged to dispose of the batteries.
- Note: Batteries containing pollutant substances are marked as follows:  
Pb= Battery contains lead,  
Cd = Battery contains cadmium,  
Hg = Battery contains mercury.
- The appliance should be disposed of according to Regulation 2002/96/EC-WEEE (Waste Electrical and Electronic Equipment). In case of queries, contact the municipal authorities responsible for waste disposal in your area.

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