
# KORONA<sup>®</sup>

## LEONIA

- D** GLASWAAGE Gebrauchsanleitung  
**GB** GLASS SCALE Instruction for use  
**FR** PÈSE-PERSONNE EN VERRE Mode d'emploi  
**ES** BÁSCULA DE VIDRIO PARA PERSONAS Instrucciones para el uso  
**IT** BILANCIA IN VETRO Istruzioni per l'uso  
**GR** ΓΥΑΛΙΝΗ -ΥΓΑΡΙΑ Οδηγίες χρήσεως  
**RUS** СТЕКЛЯННЫЕ ВЕСЫ Инструкция по применению  
**PL** WAGA SZKLANA Instrukcja obsługi  
**NL** WEEGSCHAAL MET GLAZEN WEEGPLATEAU Gebruikshandleiding  
**PT** BALANÇA DE VIDRO Instruções de utilização  
**TR** CAM TERAZI Kullanma Talimatı  
**DK** GLASVÆGT Brugsanvisning  
**SV** GLASVÅG Bruksanvisning  
**FIN** LASIVAAKA Käyttöohje  
**CZ** SKLENĚNÁ VÁHA Návod k použití  
**SL** STEKLENA TEHTNICA Navodila za uporabo  
**HU** ÜVEGMÉRLEG Használati útmutató

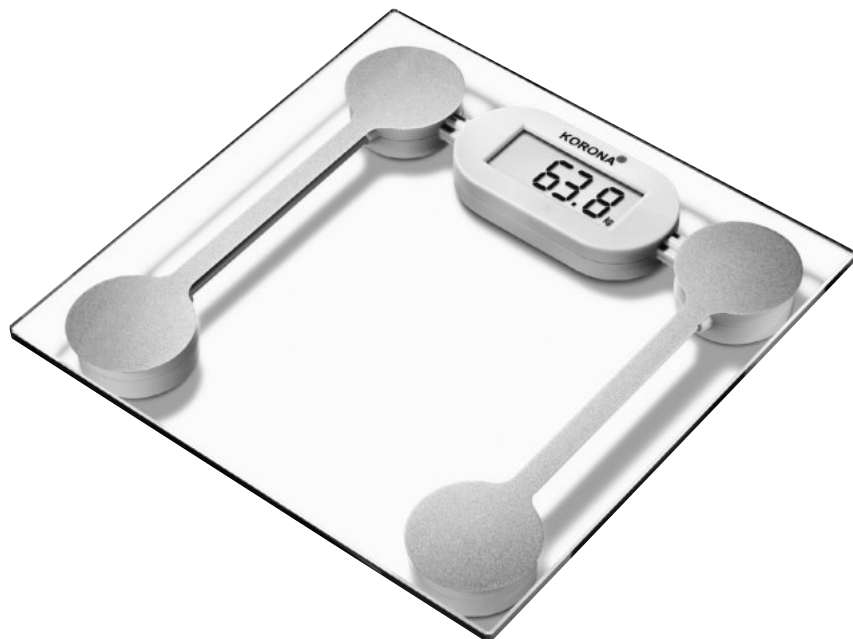Hans Dinslage GmbH  
Riedlinger Straße 28  
88524 Uttenweiler, Germany  
Tel.: +49 7374 / 915766  
[www.korona-gmbh.de](http://www.korona-gmbh.de)

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- (GB) Instructions for use  
(F) Instructions d'utilisation  
(E) Instrucciones de uso  
(I) Istruzioni per l'uso  
(GR) Οδηγίες χρήσης
- (RU) Инструкция по применению  
(PL) Instrukcja obsługi  
(NL) Gebruikershandleiding
- (P) Instruções de utilização  
(TR) Kullanım kılavuzu  
(DK) Brugsanvisning  
(S) Bruksanvisning  
(FIN) Käyttöohje
- (CZ) Návod k použití  
(SLO) Navodila za uporabo  
(H) Használati utasítás

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- (GB) Getting started  
(F) Mise en service  
(E) Puesta en marcha  
(I) Messa in funzione  
(GR) Έναρξη λειτουργίας
- (RU) Ввод в эксплуатацию  
(PL) Uruchomienie  
(NL) Ingebruikname
- (P) Colocação em funcionamento  
(TR) İlk çalıştırma  
(DK) Ibrugtagning  
(S) Börja använda vågen
- (FIN) Käyttöönotto  
(CZ) Uvedení do provozu  
(SLO) Prvi vklop  
(H) Üzembe helyezés

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### GB General advice

- **Cleaning:** You can clean the scale with a damp cloth and a little washing up liquid, if required. Never immerse the scale in water or rinse it under running water.
- Do not use aggressive or abrasive cleaning products.
- Protect the scales against knocks, damp, dust, chemicals, drastic changes in temperature, electromagnetic fields and sources of heat that are too near.
- Never step onto the outer edge of the scales with one foot as this could cause the scales to tip!
- Keep packaging material away from children!

- The device is for private use only and is not intended for medical or commercial purposes.
- Should you have any questions regarding operation and usage of our appliances, please contact your local retailers or Customer Service.
- Before making a claim, please check the batteries and replace them if necessary.
- Repairs must only be carried out by Customer Services or authorised suppliers.
- All our personal scales comply with EU guideline 2004/108/EC + relevant supplements.
- Spent batteries and rechargeable batteries do not constitute normal household waste! They are considered to be toxic waste and, as such, should be disposed of in special containers, toxic waste collection points or brought to electrical goods dealers.

You are legally obliged to dispose of the batteries.

- Note: Batteries containing pollutant substances are marked as follows:

Pb= Battery contains lead,  
Cd = Battery contains cadmium,  
Hg = Battery contains mercury.

- Follow local regulations regarding the disposal of materials. The appliance should be disposed of according to Regulation 2002/96/EC-WEEE (Waste Electrical and Electronic Equipment). In case of queries, contact the municipal authorities responsible for waste disposal in your area.

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